

HemoTemp II Heating Block Cover

The HemoTemp II activator heating block cover is an optional add-on to the HemoTemp II activator and designed to be placed on top of the heating block. This cover creates channels for the HemoTemp II indicators and prevents them from falling sideways into the block.


Align the cover on top of the HemoTemp II activator heating block.



Ensure the half circle portion of the cover is positioned around the block retrieval hole, leaving this hole unobstructed.



Note: The cover has pins protruding from the two center columns to prevent movement and create a snug fit.

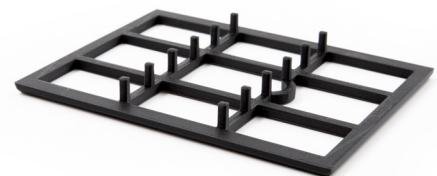

Place indicators into the heating block.

The cover does not need to be removed from the block and can stay attached during heating cycles and when the activator is turned off.

